Lifenaxx[®] * User Manual

Lifenaxx Exercise and mobility trainer LX-023

Before using the device, the first time, please read the user manual carefully. Service phone No. for technical support: 01805 013643 (14 cent/minute from German fixed-line and 42 cent/minute from mobile networks).

Free Email: support@technaxx.de

Keep this user manual for future reference or product sharing carefully. Do the same with the original accessories for this product. In case of warranty, please contact the dealer or the store where you bought this product.

Warranty 2 years

*Enjoy your product. *Share your experience and opinion on one of the well-known internet portals.*

Features:

TIME -----00: 00-99:59MIN

COUNT-----0-9999

REPS/MIN-----0-9999

CALORIE -----0-9999KAL

KEY FUNCTION:

MODE: This key lets you to select and lock on to a particular function you want. Hold the key for 4 seconds to have all function values reset

AUTO ON/OFF: The system turns on when any key is pressed or when exercise start. The monitor will be automatically shut off if there is no use after 4 minutes.

TIME By pressing the MODE key until the „>“ points to TIME. When the speed sensor has any signal input, it starts the count up and automatically accumulates workout time when starting exercise.

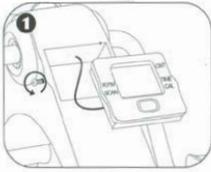
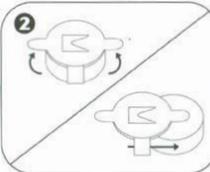
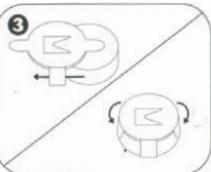
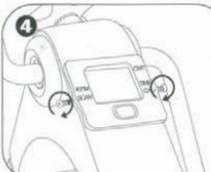
CNT show the number of rotations

RPM: show the actual speed of rotation by minute

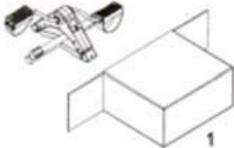
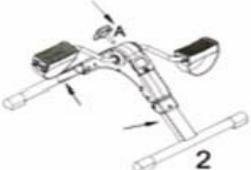
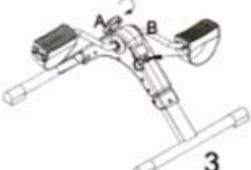
CAL: shows the actual calories burned during workout

SCAN: display's the different informations: Time, CNT, CAL, RPM every 5 seconds the display change

Change the Battery:

			
Remove the Display	Open the battery holder	Change the Battery (+ up) Typ LR044	Set in the Display

First use:

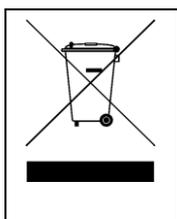
		
Unpack and remove all packing, unfold the LX-023, make sure that the spring bolts lock securely	Fix the resistance screw to your preferred resistance level	Start your exercise on arms or legs

Make sure the item stands on even Floor or desk, if the surface is too slippery place an anti-slip mat under the Exercise and mobility trainer.

Package content	Lifenaxx LX-023, user manual
-----------------	------------------------------

Warnings

- Do not stand on the exercise and mobility trainer
- The item is made to train while seated only.
- Before each workout ensure the item is unfolded correctly and the spring bolt are fixed securely.
- While using Exercise and mobility trainer make sure the underground is even and item stands stable.
- Wear appropriate footwear while performing exercise
- Perform exercises in a slow and controlled manner.
- These guidelines are here to protect you.
- We do not accept liability from injuries due to the use or misuse of this product.

Hints for Environment Protection: Packages materials are raw materials and can be recycled. Do not disposal old devices or batteries into the domestic waste.

Cleaning: Protect the device from contamination and pollution (use a clean drapery). Avoid using rough, coarse-grained materials or solvents/aggressive cleaner. Wipe the cleaned device accurately. **Distributor:** Technaxx Deutschland GmbH & Co.KG, Kruppstr. 105, 60388 Frankfurt a.M., Germany

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