

Lifenaxx Neck massager LX-046

Before using the appliance for the first time, please read the instructions for use and safety information carefully.



This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capacities, or by persons lacking in experience or knowledge, unless they are supervised or instructed on the use of this device by a person responsible for their safety. Children should be supervised to ensure they do not play with this device.

Keep this user manual for future reference or product sharing carefully. Do the same with the original accessories for this product. In case of warranty, please contact the dealer or the store where you bought this product.

Enjoy your product. * Share your experience and opinion on one of the well-known internet portals.

Specifications are subject to change without notice - please be sure to use the latest manual available on the manufacturer's website.

Table of Contents

Hints	2	Charging.....	5
Safety instructions	2	Turn device On/Off.....	5
Safety rules for rechargeable batteries.....	3	Technical Specifications.....	5
Warnings	4	Support.....	5
Disclaimer.....	4	Care and maintenance	5
Features	4	Declaration of Conformity.....	6
Product overview.....	5	Disposal.....	6

Hints

Only use the product for the purposes due to its intended function. Do not damage the product. Following cases may damage the product: Accidents (including liquid or moisture), improper or abuse use of the product, insects, tampering or modification of the product by persons other than authorized service personnel, exposure to abnormally corrosive materials, insertion of foreign objects into the unit, used with accessories not preapproved.

Read and observe all warnings and precautions in the operating instructions. Refer to and heed all warnings and precautions in the user manual. Children must not play with the device. If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of. The device has a hot surface. Persons insensitive to heat must be careful when using the device.

Safety instructions

- Never touch the product with wet or damp hands.
- The cable is sheathed and insulates the current. Do not place any objects on the cable, as this could cause the protective function to be lost. Avoid placing excessive stress on the cable, especially at sockets and plugs.
- If you want to pull the plug out of the socket, always pull on the plug and not on the cable. The cable could otherwise break.
- Do not overload wall outlets, extension cords, or other wiring as this may result in fire or electric shock.
- Do not use liquid cleaners or cleaning sprays. The product may only be cleaned with a damp cloth.
- Do not use the product in direct sunlight or places where the temperature may exceed 50°C for a long time.

●Keep packaging material away from children. There is the risk of suffocation.

Do not use the massager:

- on animals
- if you suffer from a medical abnormality or injury in the area to be massaged (e.g. herniated disc, open wounds)
- under any circumstances for massages in the area around the heart if you wear a pacemaker – and consult your doctor regarding other parts of the body
- on swollen, burnt, inflamed or damaged skin and body parts of any kind
- on cuts, capillaries, varicose veins, acne, couperose, herpes or other skin diseases
- on the head or on hard parts of the body (e.g. bones and joints)
- during pregnancy
- on the face (eyes), on the larynx or other particularly sensitive body parts
- for foot massages
- while sleeping
- after the consumption of medicines or alcohol (reduced cognition!)
- on skin with previously applied cream (e.g. warming ointment),
- in vehicles.

Consult your doctor before using the massager, especially:

- if you suffer from a severe illness or have recently had surgery on your upper body
- if you have a pacemaker, implants or other aids
- if you suffer from thrombosis
- if you suffer from diabetes
- if you are suffering pain from unknown causes

Safety rules for rechargeable batteries

The product uses rechargeable batteries and can go wherever you need them. But ALL batteries can EXPLODE, CATCH FIRE, and CAUSE BURNS if disassembled, punctured, cut, crushed, short-circuited, incinerated, or exposed to water, fire, or high temperatures, so you must handle them with care.

To use the rechargeable batteries safely, follow these guidelines:

- NEVER let the rechargeable battery or battery compartment get wet.

● NEVER disassemble, cut, crush, puncture, short-circuit, dispose of batteries in fire or water, or expose the rechargeable battery to temperatures above 50°C.

Warnings

● Do not disassemble the product, it may result in short-circuit or damage.
● Battery life-span will decrease with usage. Charge the battery at least once every month to preserve battery in the event the product is not used for long time. ● Battery warning: Improper use of the battery may cause fire or chemical burns. The battery may explode in case of damage. ● Do not modify, repair or remove without professional guidance. ● Do not use corrosive or volatile liquid for cleaning. ● Do not drop or shake the device, it may break internal circuit boards or mechanics. ● Keep the device away from little children.

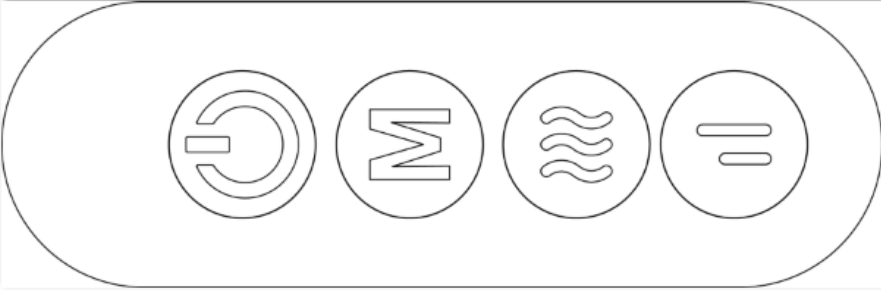




Disclaimer

● In no event shall Technaxx Deutschland be liable for any direct, indirect punitive, incidental, special consequential danger, to property or life, improper storage, whatsoever arising out of or connected with the use or misuse of their products.

Features

● The massager adapts easily to different parts of the body.
● The soothing heat can be switched off.
● You determine the intensity and course of the massage yourself.
● With 2 different modes and 2 different adjustable speeds, it increases well-being, relieves pain, sore muscles, cramps and muscle tension.
● The massager is suitable for use on the neck, shoulders, back, waist, abdomen, arms and legs.
● 2 modes
● 2 different speed

Product overview

Control panel			
			
	On/ Off switch		Program
	Heat (LED red ->on, LED green -> off)		Speed

Charging

●Connect the USB-C cable on the device to a USB adapter* (*not included).

Note: First charging takes about 3h, RED LED flash as long as item load.

Turn device On/Off

Position the massager on a chosen place on your body

Press and hold the power button until it lights up red

Technical Specifications

Battery capacity	2000mAh 7.4V
Rated power	10W
Wireless use, with USB cable and Type-C charging port	No use during charging
Automatic shut-off function after 10 minutes	
Heating inside could be switched off.	

Support

Service phone No. for technical support: **01805 012643** (14 cent/minute from German fixed-line and 42 cent/minute from mobile networks). Free

Email: **support@technaxx.de**

The support hotline is available Mon-Fri from 9am to 1pm & 2pm to 5pm

Care and maintenance

Clean the device only with a dry or slightly damp, lint-free cloth.

Do not use abrasive cleaners to clean the device.

- Use the device in ultra-high or ultra-low temperature.
- Deliver or use it in strongly shocking environment.

Declaration of Conformity



The EU Declaration of Conformity can be requested at the following address: www.technaxx.de/ (in the lower bar "Declaration of Conformity").

Disposal



Disposal of the packaging. Sort packaging materials by type upon disposal.

Dispose of cardboard and paperboard in the waste paper. Foils should be submitted for recyclables collection.



Disposing of old equipment (Applies in the European Union and other European countries with separate collection (collection of recyclable materials) Old equipment must not be disposed of with household waste! Every consumer is required by law to



dispose of old devices that can no longer be used separately from household waste, e.g. at a collection point in his or her municipality or district. This ensures that the old devices are properly recycled and that negative effects on the environment are avoided. For this reason, electrical devices are marked with the symbol shown here.



Batteries and rechargeable batteries must not be disposed of in household waste! As a consumer, you are required by law to dispose of all batteries and rechargeable batteries, whether they contain harmful substances* or not, at a collection point in your community/city or with a retailer, to ensure that the batteries can be disposed of in an environmentally friendly manner. * marked with: Cd = cadmium, Hg = mercury, Pb = lead. Return your product to your collection point with the fully discharged battery installed inside!

Made in China

Distributed by:
Technaxx Deutschland GmbH & Co. KG
Konrad-Zuse-Ring 16-18,
61137 Schöneck, Germany

Neck massager LX-046