

Lifenaxx[®] * user manual

Micro Massager LX-039

Before using the appliance for the first time, please read the instructions for use and safety information carefully.



Keep this user manual for future reference or product sharing carefully. Do the same with the original accessories for this product. In case of warranty, please contact the dealer or the store where you bought this product.

Enjoy your product * Share your experience and opinion on one of the well-known internet portals.

Specifications are subject to change without notice - please be sure to use the latest manual available on the manufacturer's website.

Features

The Lifenaxx LX-039 is a cordless vibration massager that relaxes sore, stiff muscles, increases circulation and range of motion, and improves the overall health of the body's soft tissues.

- Speed up to 3200 rpm.
- 7.4V Volt battery with 2000mAh
- Charging time approx. 2 hours
- Operating time approx. 2 hours
- 4 different massage heads
- Battery indicator


Technical specifications

USB connection	USB-C with 5V, max. 3A
Battery	7.4V, 2000mAh
Power	25W
Speed	Max. 3200rpm
Duration	Approx. 2 hours
Loading time	Approx. 2 hours
Engine type	Brush motor
Material	ABS plastic
Accessories	LX-039 Micro Massager, 4x message heads, charging cable, user manual

Usage

● Charge indicator light: The blue indicator light flashes slowly when charging. Fully charged when all lights light up. ● Install 1 of the 4 heads according to your needs. ● Press the power button to switch the unit on/off. ● Press the button briefly to change the speed. ● Install one of the 4 heads according to your needs. ● Use only over clothing, on dry, clean body surface, by moving the device with light pressure for about 60 seconds per body area. ● Use the Lifenaxx LX-039 only on the soft tissue of the body, as needed, without causing pain or discomfort. Do not use on the head or on hard or bony areas of the body. ● Only use the applicators that give the best possible effect. ● Hematomas may occur, regardless of the pressure applied. Check the treated areas frequently and stop immediately if you notice the first signs of pain or discomfort. ● Keep fingers, hair or other body parts away from the shaft and back of the applicator as this may cause bruising.

Product overview

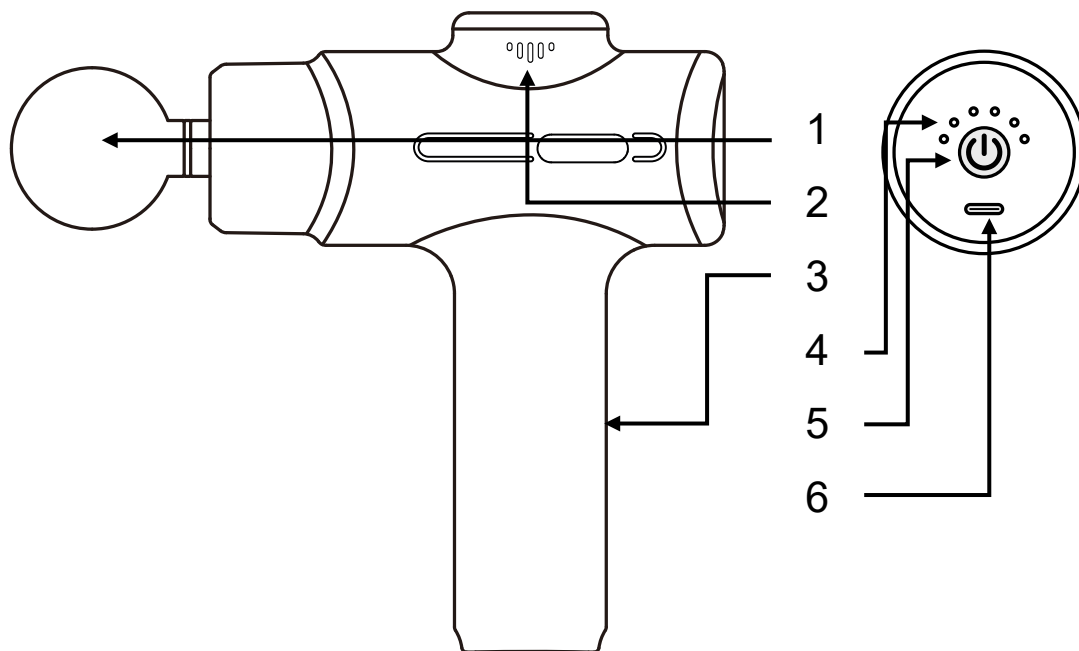

1	Massage head	4	Battery LED indicator
2	Air outlet	5	Start/Stop button
3	Handle	6	USB Type-C charging port

Attachments from left to right



Fork head	Suitable for the spine
Pin head	Suitable for parts close to the joints and deep stimulation
Flat head	Suitable for large muscle areas
Soft head	Suitable for muscle groups

Installing the attachments

Install and replace the massage head as needed. To do this, simply pull off the attachment and insert the desired attachment into the tapered round hole of the massager.

Safety instructions

- Only charge the unit with the USB charging cable supplied.
- This appliance is not intended for children or persons with reduced mental capacity or lack of experience and/or expertise.
- This appliance is not intended for commercial use.
- Any other use or modification of the appliance is considered improper and involves considerable risks. The manufacturer accepts no liability for damage caused by improper use.
- Never leave the Lifenaxx LX-039 unattended when in operating mode or when charging.

PLEASE DO NOT USE THE LIFENAXX LX-039 IN THE FOLLOWING CIRCUMSTANCES WITHOUT THE PRIOR CONSENT OF YOUR DOCTOR:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wearing pacemakers, recent surgery, epilepsy or migraine, herniated disc, spondylolisthesis, spondylolysis or spondylosis, recent joint replacement surgery or IUDs, metal pins or plates or if you have any concerns about your physical health. These contraindications do not mean that you cannot use a vibration or exercise machine. However, we recommend that you consult a doctor first. The effects of vibration use on certain medical conditions are currently being investigated. This will probably lead to a shortening of the list of contraindications given. Practice has shown that the integration of vibration applications into a treatment plan is advisable in a number of cases. This must be done on the recommendation and with the guidance of a doctor, specialist or physiotherapist.

Warnings & Caution

- Do not use the device with life support equipment if failure of the device may result in failure of the life support equipment or significant reduction in safety or effectiveness.
- Do not use the device in the presence of a flammable anesthetic mixture with air, oxygen or nitrous oxide, or near flammable substances, vapors or gases.
- Do not operate the unit near heating vents or in direct sunlight.
- Keep the unit dry and switch it off when not in use.

Support

Service phone No. for technical support: **01805 012643*** (14 cent/minute from German fixed-line and 42 cent/minute from mobile networks). Free

Email: **support@technaxx.de**

The support hotline is available Mon-Fri from 9am to 1pm & 2pm to 5pm.

Declaration of Conformity



The EU Declaration of Conformity can be requested at the following address: www.technaxx.de/ (in the lower bar "Declaration of Conformity").

Disposal



Disposal of the packaging. Sort packaging materials by type upon disposal.

Dispose of cardboard and paperboard in the waste paper. Foils should be submitted for recyclables collection.



Disposing of old equipment (Applies in the European Union and other European countries with separate collection (collection of recyclable materials) Old equipment must not be disposed of with household waste! Every consumer is required by law to dispose of old devices that can no longer be used separately from household waste, e.g. at a collection point in his or her municipality or district. This ensures that the old devices are properly recycled and that negative effects on the environment are avoided. For this reason, electrical devices are marked with the symbol shown here.



Batteries and rechargeable batteries must not be disposed of in household waste! As a consumer, you are required by law to dispose of all batteries and rechargeable batteries, whether they contain harmful substances* or not, at a collection point in your community/city or with a retailer, to ensure that the batteries can be disposed of in an environmentally friendly manner. * marked with: Cd = cadmium, Hg = mercury, Pb = lead. Return your product to your collection point with the fully discharged battery installed inside!

Made in China

Distributed by:

Technaxx Deutschland GmbH & Co. KG
Konrad-Zuse-Ring 16-18,
61137 Schöneck, Germany